




# **CLICK-&-GO® PLUS 555 OPP PROTEIN SYNTHESIS ASSAY KIT**

**SKU:** CCT-1494

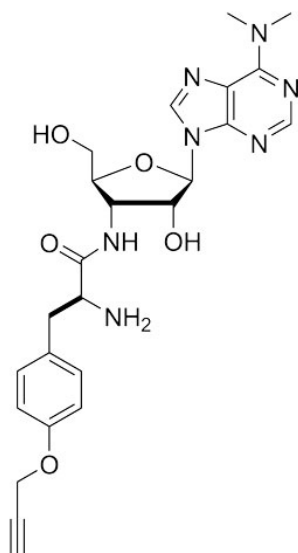

---

## **DESCRIPTION**

Synthesis of many proteins is tightly controlled at the level of translation, and plays an essential role in fundamental processes such as cell growth and proliferation, signaling, differentiation, or death. Methods that allow imaging and identification of nascent proteins are critical for dissecting regulation of translation, both spatially and temporally, particularly in whole organisms. Although protein synthesis is a conserved and essential cellular function, it is often regulated in a cell-type-specific manner to influence cell fate, growth and homeostasis. Most methods used to measure protein synthesis depend on metabolically labeling large numbers of cells with radiolabeled amino acids, stable isotope-labeled amino acids, bioorthogonal noncanonical amino acid tagging (L-azidohomoalanine or homopropargylglycine or their combination). Because these methods typically depend on specialized growth conditions, they

**For research use only. Not intended for therapeutic or diagnostic use in animals or humans.**



have been largely restricted to yeast, bacteria and cell lines. Application of these techniques for investigating protein synthesis within mammalian systems *in vivo* has been challenging.


The use of O-propargyl-puromycin (OPP), an analog of puromycin that contains a terminal alkyne group, has facilitated the quantification of protein synthesis within individual cells *in vivo*. OPP enters the acceptor site of ribosomes and incorporates into nascent polypeptide chains. Unlike traditional methods mentioned above, OPP is not an amino acid analog; thus, OPP can be added directly to cells in complete media (i.e., methionine-containing) or used to detect *in vivo* protein synthesis. It also can be used with cell lines that are sensitive to media exchanges or incubation in methionine-free media. The combination of high cell permeability and signal-to-noise ratio makes OPP an ideal candidate compound to study nascent proteomes across a wide array of cellular types and conditions. The kit contains all of the components needed to detect incorporated OPP with orange-fluorescent AZDye 555 Azide Plus (Alexa Fluor® 555 equivalent), and blue-fluorescent Hoechst 33342 dye for nuclear staining. A sufficient amount of reagents is provided for imaging 25 coverslips or 250 wells using 96-well plates.

## SPECIFICATIONS

<b>Unit Size</b>	3 Kit
<b>Storage Instructions</b>	4C
<b>Excitation/Emission Maximum</b>	556/572 nm
<b>Label</b>	AZDye 555 Azide Plus
<b>Number of Reactions</b>	25
<b>Shipping Instructions</b>	Ambient temperature

## ABS/EM SPECTRA


**For research use only. Not intended for therapeutic or diagnostic use in animals or humans.**

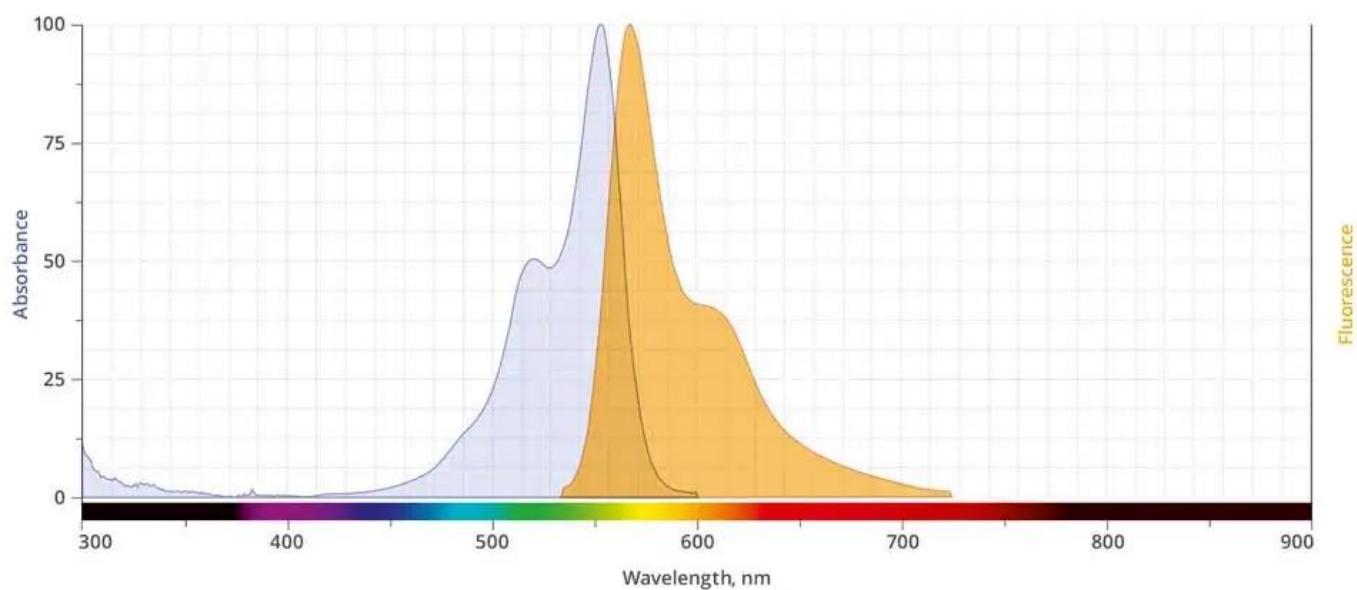

## DOCUMENTS

- [Safety Data Sheet](#)
- [Download CoA](#)
- [User Guide](#)
- [Datasheet](#)

## GALLERY IMAGES

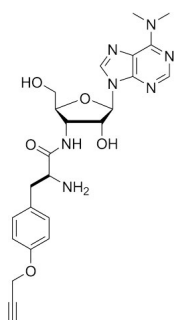

**For research use only. Not intended for therapeutic or diagnostic use in animals or humans.**