




5-TAMRA AZIDE

SKU: CCT-1245



DESCRIPTION

5-TAMRA Azide (also known as TAMRA azide, isomer 5) is the red-fluorescent probe that is compatible with various excitation sources including mercury arc, tungsten and xenon arc lamps, the 544 nm line of the Helium-Neon laser and the 532 nm green laser line. It is predominantly used for detection of terminal alkyne-tagged biomolecules via a copper-catalyzed click reaction (CuAAC). It also reacts with strained cyclooctyne via a copper-free “click chemistry” reaction to form a stable triazole and does not require Cu-catalyst or elevated temperatures.

Although the mixed isomers of 5(6)-TAMRA Acid is a preferred, routinely used orange-fluorescent dye for staining proteins, it is also used for labeling peptides and nucleotides. Purification of 5(6)-TAMRA labeled peptide and nucleotides might be troublesome due to significant signal broadening in HPLC purification. Peptides and nucleotides labeled with a single isomer TAMRA usually give better resolution in HPLC purification that is often required in the conjugation processes. 5-TAMRA is more often used for labeling peptides and proteins. 6-TAMRA is predominately used for labeling nucleotides and sequencing nucleic acids. 5-TAMRA dye is a bright fluorescent label is compatible with various excitation sources including mercury arc,

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.



tungsten and xenon arc lamps, the 544 nm line of the Helium-Neon laser and the 532 nm green laser line.


TAMRA NHS Ester reacts specifically and efficiently with a primary amine (e.g., side chain of lysine residues or aminosilane-coated surfaces) at pH 7-9 to form a stable, covalent amide bond. The NHS ester (or succinimidyl ester) is the most popular tool for conjugating dyes to the primary amines of protein or antibody (Lys), amine-modified oligonucleotides, and other amine-containing molecules.

SPECIFICATIONS

CAS Number	150810-68-7
Molecular Weight	527.53
Appearance	Dark red amorphous solid
Extinction Coefficient	92,000
Purity	>90% (HPLC)
Unit Size	1 mg, 5 mg, 25 mg, 100 mg
Solubility	Water, DMSO, DMF
Storage Instructions	-20°C. Desiccate
Spectrally Similar Dyes	Alexa Fluor® 555, CF® 555, DyLight® 549
Laser Line	532 nm or 555 nm
Excitation/Emission Maximum	553/575 nm
Shipping Conditions	Ambient temperature
Shipping Instructions	Ambient temperature

ABS/EM SPECTRA


For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

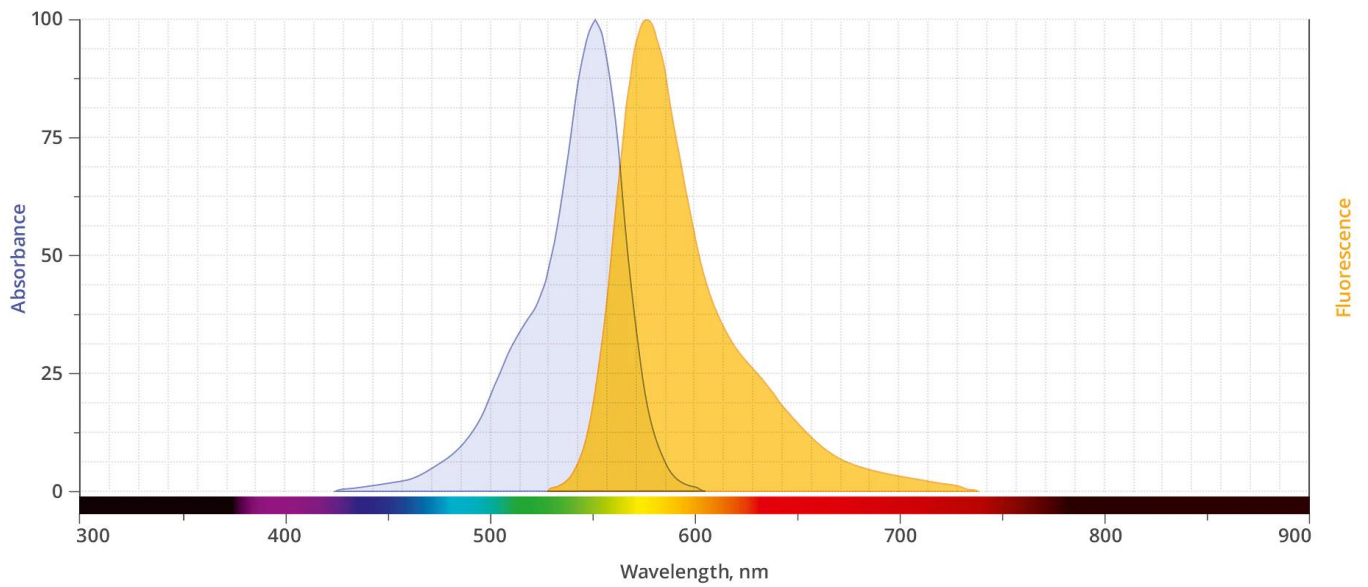

DOCUMENTS

- [Safety Data Sheet](#)
- [Download CoA](#)
- [Datasheet](#)

GALLERY IMAGES



For research use only. Not intended for therapeutic or diagnostic use in animals or humans.