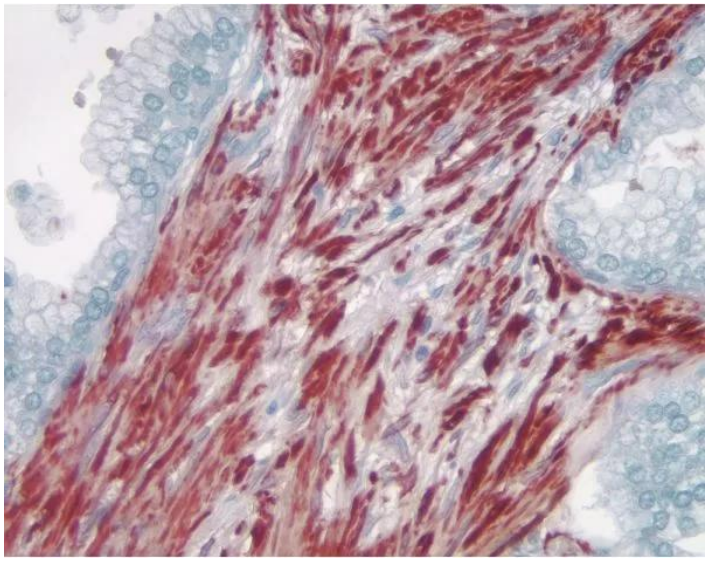
www.vectorlabs.com
Email: customerservice@vectorlabs.com
Telephone: (650) 697-3600

AEC SUBSTRATE KIT, PEROXIDASE (HRP), (3-AMINO-9-ETHYLCARBAZOLE)

SKU: SK-4200



DESCRIPTION

AEC (3-amino-9-ethylcarbazole) HRP substrate produces a red reaction product that can be used as a single label or as a second color for multiple antigen labeling. Stained slides must be aqueously mounted because these reaction products are soluble in organic solvents.

Features:

- Greater sensitivity than conventional substrates
- Consistent and reliable
- Suitable for IHC/ICC and blots
- Stock solutions supplied in convenient dropper bottles promoting ease of handling
- No wait times for mixing and dissolving powders or tablets
- Sufficient reagents to produce 300 ml of working solution

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

www.vectorlabs.com
Email: customerservice@vectorlabs.com
Telephone: (650) 697-3600

SPECIFICATIONS

Chromogen Color	Red
Detection Enzyme	Peroxidase
Microscopy	Brightfield, Spectral Imaging
Mounting Medium Type	Aqueous (Hardening), Aqueous (Non-Hardening)
Unit Size	1 Kit
Applications	Immunohistochemistry / Immunocytochemistry, In situ hybridization, Blotting Applications, Elispot
Additional Substrate Properties	Contrast with Tissue Pigments, Multiple Labeling
Safety Title	WARNING: Cancer –www.P65Warnings.ca.gov
Substrate Reaction Product	Precipitating Reaction Product

KIT COMPONENTS

- 6 ml AEC Reagent 1
- 6 ml AEC Reagent 2
- 6 ml AEC Reagent 3

TECHNICAL INFORMATION

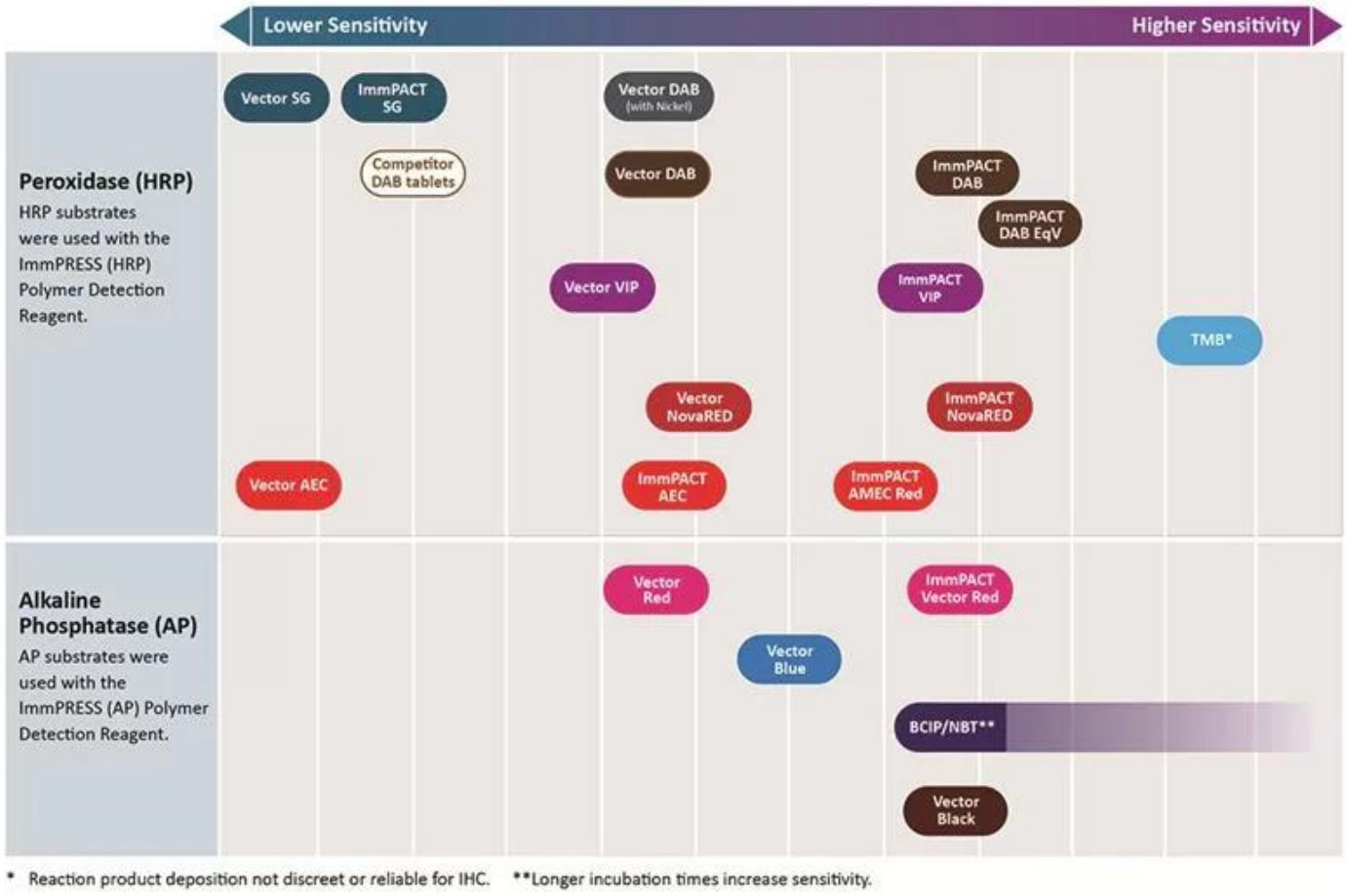
AEC and ImmPACT AEC chromogens can be used for both manual and automated staining methods. Sections stained with AEC are stable for at least 4 years when mounted in VectaMount® AQ (H-5501).

This kit contains stock solutions in convenient dropper bottles. This product can also be used on blots. The kit provides all of the necessary reagents to prepare about 300 ml of working solution. When stored at 4 C, this kit is stable for one year.

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

www.vectorlabs.com
Email: customerservice@vectorlabs.com
Telephone: (650) 697-3600

Relative Sensitivity of Substrates in Immunohistochemistry



* Reaction product deposition not discreet or reliable for IHC. **Longer incubation times increase sensitivity.

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

www.vectorlabs.com
Email: customerservice@vectorlabs.com
Telephone: (650) 697-3600

Counterstain/Substrate Compatibility Table

This table is designed as a reference to determine the optimal counterstain/substrate combination for your application. Considerations should be given to tissue type, antigen unmasking protocol and other detection parameters to achieve the desired staining intensity.

Substrate	Catalog Number	Vector Hematoxylin and Hematoxylin QS H-3401 and H-3404	Vector Methyl Green H-3402	Vector Nuclear Fast Red H-3403
ImmPACT DAB (brown)	SK-4105	Excellent Contrast	Excellent Contrast	Fair Contrast
ImmPACT DAB EqV	SK-4103	Excellent Contrast	Excellent Contrast	Fair Contrast
DAB (brown)	SK-4100	Excellent Contrast	Excellent Contrast	Fair Contrast
DAB-Ni (gray-black)	SK-4100	Excellent Contrast	Fair Contrast *	Good Contrast
ImmPACT AEC (red)	SK-4205	Excellent Contrast	Counterstain Incompatibility **	Color Incompatibility
ImmPACT AMEC Red (red)	SK-4285	Excellent Contrast	Counterstain Incompatibility **	Color Incompatibility
AEC (red)	SK-4200	Excellent Contrast	Counterstain Incompatibility **	Color Incompatibility
TMB (blue)	SK-4400	Color Incompatibility	Counterstain Incompatibility	Excellent Contrast
ImmPACT VIP (purple)	SK-4605	Fair Contrast	Excellent Contrast	Poor Contrast
Vector VIP (purple)	SK-4600	Fair Contrast	Excellent Contrast	Poor Contrast
ImmPACT SG (blue-gray)	SK-4705	Poor Contrast	Good Contrast	Excellent Contrast
SG (blue-gray)	SK-4700	Poor Contrast	Good Contrast	Excellent Contrast
ImmPACT NovaRED (red)	SK-4805	Excellent Contrast	Excellent Contrast ***	Color Incompatibility
Vector NovaRED (red)	SK-4800	Excellent Contrast	Excellent Contrast ***	Color Incompatibility
ImmPACT Vector Red (magenta)	SK-5105	Excellent Contrast	Excellent Contrast	Color Incompatibility
Vector Red (magenta)	SK-5100	Excellent Contrast	Excellent Contrast	Color Incompatibility
Vector Black (black)	SK-5200	Excellent Contrast	Excellent Contrast *	Excellent Contrast
Vector Blue (blue)	SK-5300	Color Incompatibility	Good Contrast	Excellent Contrast
BCIP/NBT (indigo)	SK-5400	Color Incompatibility	Excellent Contrast *	Excellent Contrast

* This substrate shows a slight decrease in sensitivity following the methyl green protocol. This decrease can be minimized by reducing the heat incubation and acetone rinse times in the methyl green protocol.

** Substrate dissolves in acetone wash.

*** A slight color change in ImmPACT NovaRED™ and Vector NovaRED™ reaction product may be seen using methyl green.

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

www.vectorlabs.com
Email: customerservice@vectorlabs.com
Telephone: (650) 697-3600

Enzyme Substrate Combinations

Recommended combinations of substrates and the recommended order in which they should be used.

First Substrate \ Second Substrate	Alkaline Phosphatase			Peroxidase					
	ImmPACT Vector Red & Vector Red (magenta) SK-5105 & SK-5100	Vector Blue (blue) SK-5300	BCIP/NBT (indigo) SK-5400	ImmPACT VIP & Vector VIP (purple) SK-4605 & SK-4600	ImmPACT DAB, ImmPACT DAB EqV & DAB (brown) SK-4105, SK-4103, SK-4100	DAB-Ni (gray-black) SK-4100	ImmPACT NovaRED & Vector NovaRED (red) SK-4805 & SK-4800	ImmPACT SG & SG (blue-gray) SK-4705 & SK-4700	ImmPACT AEC, ImmPACT AMEC Red & AEC (red) SK-4205, SK-4285, SK-4200
Alkaline Phosphatase									
ImmPACT Vector Red & Vector Red (magenta) SK-5105 & SK-5100		–	–	–	+	+	–	+	–
Vector Blue (blue) SK-5300	+		–	+	+	+	+	+	+
BCIP/NBT (indigo) SK-5400	+	–		+	+	+	+	+	+
Peroxidase									
ImmPACT VIP & Vector VIP (purple) SK-4605, SK-4600	–	+	–		+	+	–	+	–
ImmPACT DAB, ImmPACT DAB EqV & DAB (brown) SK-4105, SK-4103, SK-4100	+	+	+	+		–	–	+	+
DAB-Ni (gray-black) SK-4100	+	–	–	+	+		+	–	–
ImmPACT NovaRED & Vector NovaRED (red) SK-4805, SK-4800	–	+	+	–	+	+		+	–
ImmPACT SG & SG (blue-gray) SK-4705, SK-4700	+	–	–	+	+	–	–		+
ImmPACT AEC, ImmPACT AMEC Red & AEC (red) SK-4205, SK-4285, SK-4200	–	–	–	–	+	–	–	+	

+ Indicates good contrast – Indicates incompatiblity of substrates for various reasons

CITATIONS

Powered by Bioz © 2023 See more details on Bioz

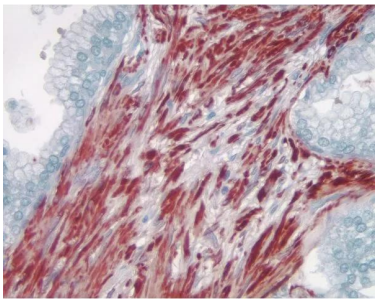
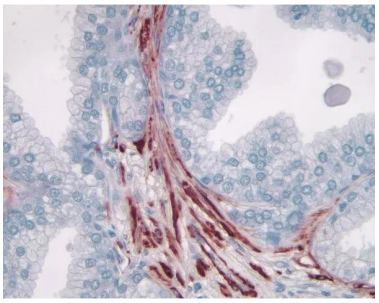
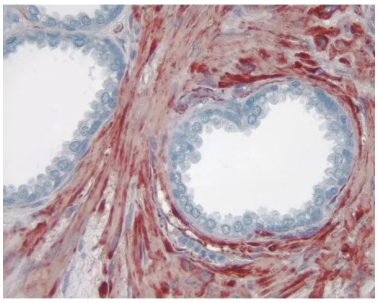
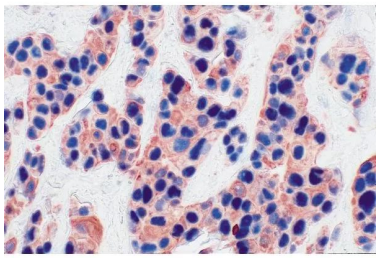
For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

www.vectorlabs.com
Email: customerservice@vectorlabs.com
Telephone: (650) 697-3600

DOCUMENTS

- Safety Data Sheet
- Enzyme Substrates Brochure
- Enzyme Substrate Combinations Chart
- Counterstain/Substrate Compatibility Chart
- User Guide
- Download CoA
- Datasheet

GALLERY IMAGES


For research use only. Not intended for therapeutic or diagnostic use in animals or humans.