




AZDYE 633 PICOLYL AZIDE

SKU: CCT-1549

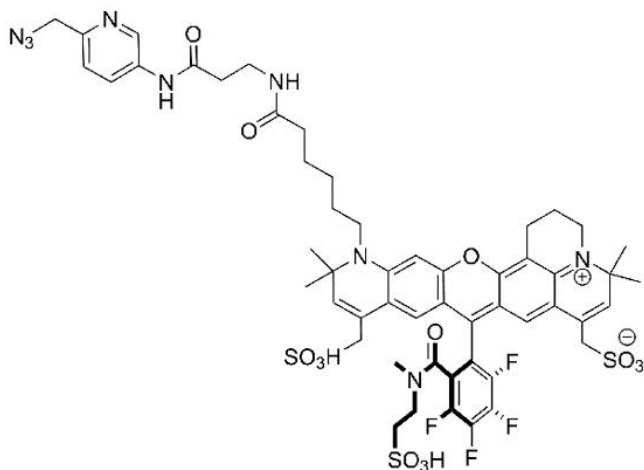

DESCRIPTION

AZDye™ 633 Picolyl Azide is an advanced fluorescent probe that incorporates a copper-chelating motif to raise the effective concentration of Cu(I) at the reaction site to boost the efficiency of the CuAAC reaction, resulting in a faster and more biocompatible CuAAC labeling. Up to 40-fold increase of signal intensity, compared to conventional azides, was reported (see Selected References).

AZDye™ 633 Azide is a bright and photostable far-red fluorescent probe with excitation ideally suited to the 633 nm or 635 nm laser excitation source. AZDye™ 633 Azide is water-soluble, pH-insensitive from pH 4 to pH 10. AZDye™ 633 spectrally is almost identical to Alexa Fluor® 633, DyLight® 633 or CF® 633 Dye. Combination of superior brightness and very low autofluorescence background signal in most biological samples in far-red spectral region allows for very sensitive detection of alkyne-labeled biomolecules.


For research use only. Not intended for therapeutic or diagnostic use in animals or humans.



SPECIFICATIONS

CAS Number	N/A
Molecular Weight	1200.26
Appearance	Blue solid
Extinction Coefficient	100,000
Unit Size	1 mg, 5 mg, 25 mg
Solubility	Water, DMSO, DMF
Storage Instructions	-20°C. Desiccate
Spectrally Similar Dyes	Alexa Fluor® 633, CF® 633
Laser Line	633 nm or 647 nm
Excitation/Emission Maximum	631/651 nm
Shipping Conditions	Ambient temperature
Shipping Instructions	Ambient temperature

ABS/EM SPECTRA

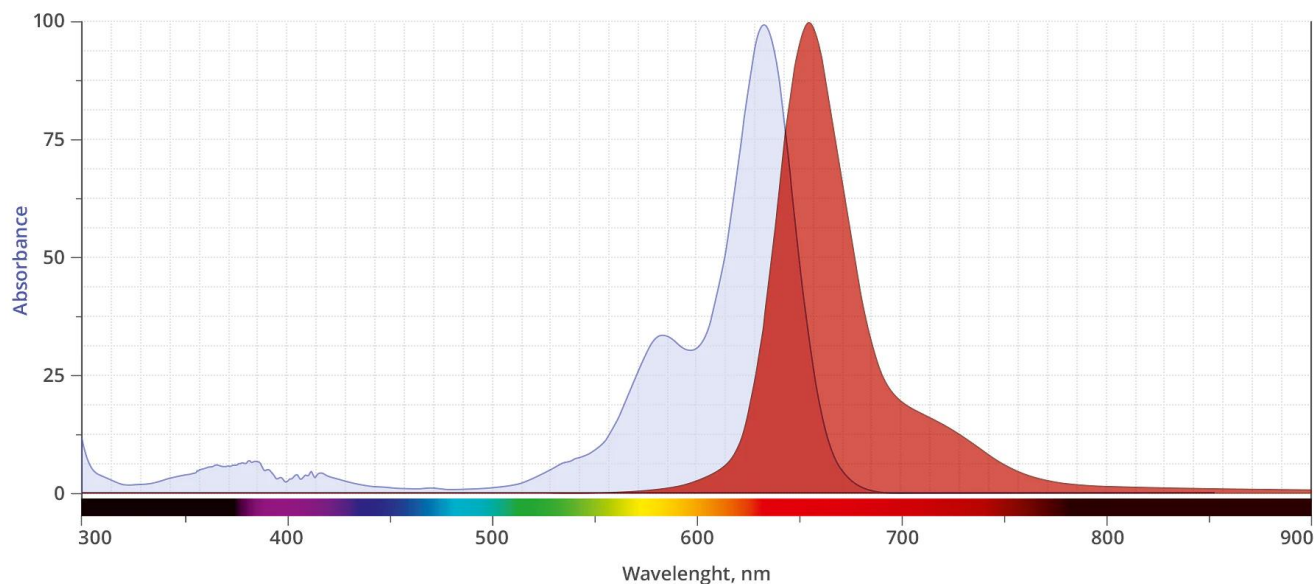

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.




SELECTED REFERENCES

1. Morral, C., et al. (2020). Protocol for Efficient Protein Synthesis Detection by Click Chemistry in Colorectal Cancer Patient-Derived Organoids Grown In Vitro. *STAR Protocols*, **Volume 1**, 2 [[ScienceDirect](#)]
2. Uchiyama, J., et al. (2020). Quantitative nascent proteome profiling by dual pulse labeling with O-propargyl-puromycin and stable isotope labeled amino acids. *The Journal of Biochemistry*, **10**, 1093. [[Oxford Academic](#)]
3. Jiang, H., et al. (2014). Monitoring Dynamic Glycosylation in Vivo Using Supersensitive Click Chemistry. *Bioconjugate Chem.*, **25**, 698-706. [[PubMed](#)]
4. Uttamapinant, C., et al. (2012). Fast, Cell-Compatible Click Chemistry with Copper-Chelating Azides for Biomolecular Labeling. *Angew. Chem. Int. Ed.*, **51**, 5852-56. [[PubMed](#)]
5. Gaebler, A., et al. (2016). A highly sensitive protocol for microscopy of alkyne lipids and fluorescently tagged or immunostained proteins. *J. Lipid. Res.*, **57**, 1934-47. [[PubMed](#)]

DOCUMENTS

- [Safety Data Sheet](#)
- [Download CoA](#)
- [Datasheet](#)

GALLERY IMAGES

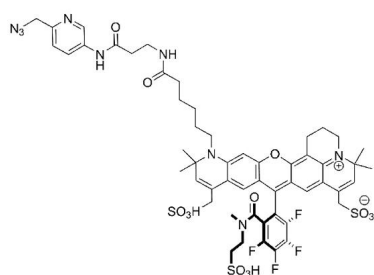

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.