




AZDYE 680 NHS ESTER

SKU: FP-1791

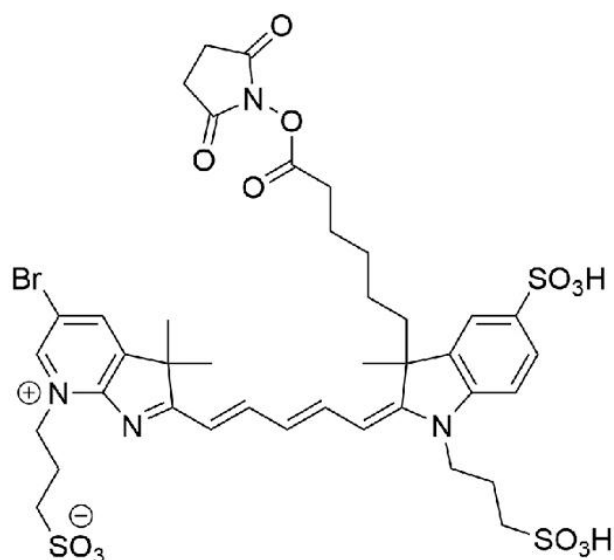

DESCRIPTION

633/647



Laser
line

Cy5.5



Common
filter set

673



Excitation
max

694



Emission
max

AZDye 680 is a bright and photostable near-IR dye that is spectrally and structurally identical to Alexa Fluor® 680 Dye. AZDye 680 is often used for small animal imaging applications in the 700

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.



nm wavelength and provides excellent clearance profiles in animals. Flow cytometry is another common application for this dye. AZDye 680 ideally suited for the 633 nm laser line. This near-IR fluorescent dye is water soluble and pH-insensitive from pH 4 to pH 10. AZDye 680 molecules can be attached to proteins at high molar ratios without significant self-quenching, leading to brighter conjugates and more sensitive detection. The long wavelength emission allows for detection in complex samples with auto-fluorescent background signals.


The NHS ester (succinimidyl ester) of AZDye 680 is the most popular tool for conjugating this dye to a protein or antibody. NHS esters can be used to label to the primary amines (R-NH₂) of proteins, amine-modified oligonucleotides, and other amine-containing molecules.

SPECIFICATIONS

Molecular Weight	955.9
Extinction Coefficient	185,000 cm ⁻¹ M ⁻¹
Reactivity	Primary amine
Unit Size	1 mg, 5 mg, 25 mg, 100 mg
Solubility	Water, DMSO, DMF
Storage Instructions	-20°C.
Spectrally Similar Dyes	Alexa Fluor® 680, DyLight® 680, IREDye® 680
Excitation/Emission Maximum	678/701 nm
Shipping Conditions	Ambient temperature
Shipping Instructions	Ambient temperature

ABS/EM SPECTRA


For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

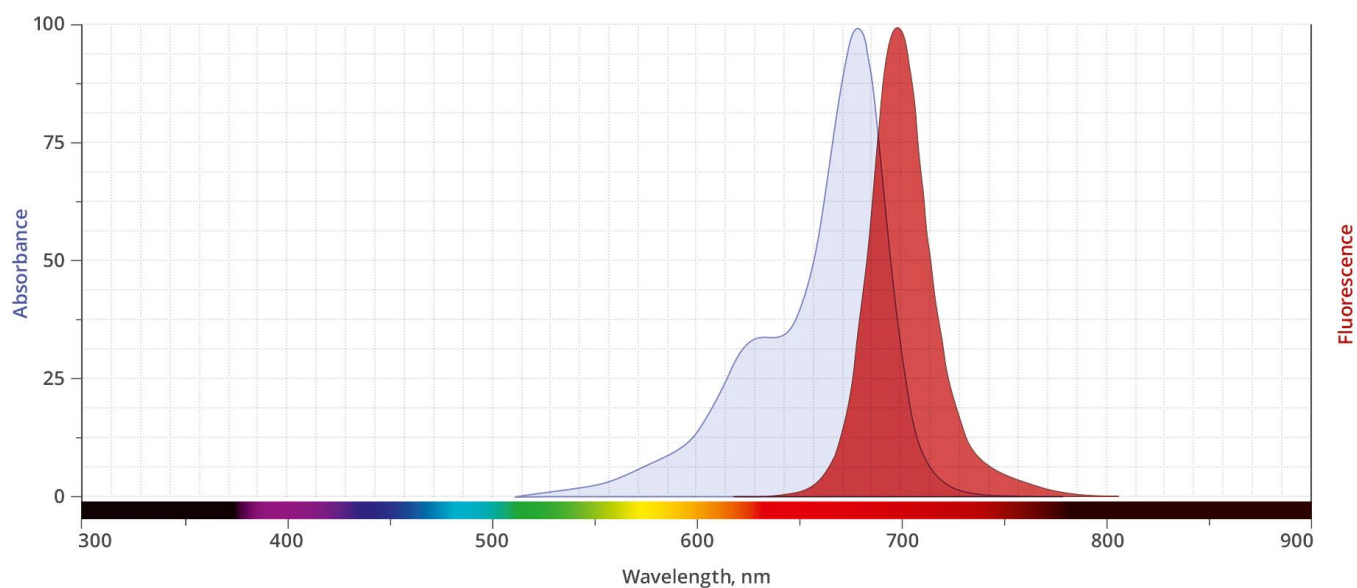

DOCUMENTS

- [Safety Data Sheet](#)
- [Datasheet](#)

GALLERY IMAGES

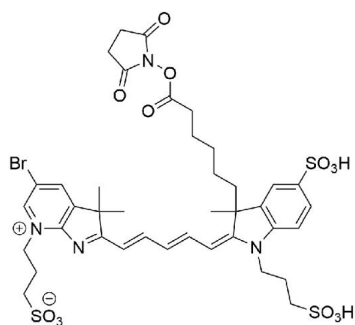

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.