




CY7 PICOLYL AZIDE

SKU: CCT-1183

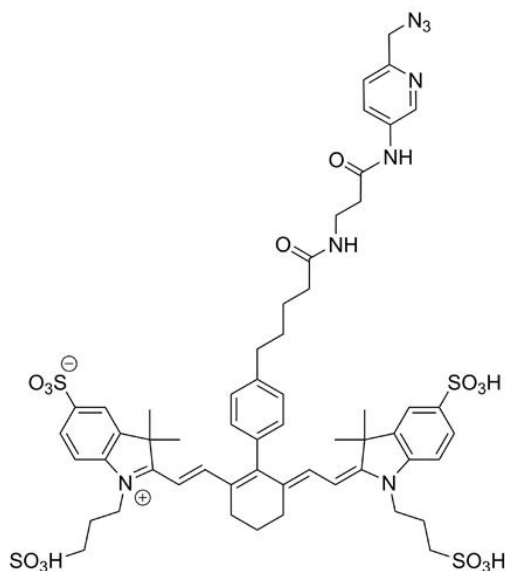

DESCRIPTION

Cy7 Picolyl Azide is an advanced fluorescent probe that incorporates a copper-chelating motif to raise the effective concentration of Cu(I) at the reaction site to boost the efficiency of the CuAAC reaction, resulting in a faster and more biocompatible CuAAC labeling. Up to 40-fold increase of signal intensity, compared to conventional azides, was reported (see Selected References).

In addition, the use picolyl azides instead of conventional azides allows for at least a tenfold reduction in the concentration of the copper catalyst without sacrificing the efficiency of labeling, significantly improving biocompatibility of CuAAC labeling protocol.

In summary, the introduction of a copper-chelating motif into azide probe leads to a substantial increase in the sensitivity and reduced cell toxicity of CuAAC detection alkyne-tagged

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.



biomolecules. This will be of special value for the detection of low abundance targets or living system imaging.


Cy7 Picolyl Azide is a bright and photostable near-IR probe that is spectrally similar to Alexa Fluor® 750, DyLight® 750, and IRDye® 750 dye. Cy7 Picolyl Azide is water-soluble, hydrophilic dye often a reagent of choice for assay where minimal non-specific binding and exceptional brightness is required. The fluorescence of Cy7 dye is pH insensitive from pH 4 to pH 10 and produces minimal autofluorescence of biological specimens in this region of the spectrum. Fluorescence of this long-wavelength Cyanine dye is not visible to the human eye but is readily detected by most imaging systems.

SPECIFICATIONS

CAS Number	N/A
Molecular Weight	1507.01
Molecular Formula	C74H111N11O14S4
Appearance	Dark green solid
Chemical Formula	C74H111N11O14S4
Extinction Coefficient	255,000
Purity	>95% (HPLC)
Unit Size	1 mg, 5 mg, 25 mg
Solubility	Water, DMSO, DMF
Storage Instructions	-20°C. Desiccate
Spectrally Similar Dyes	Alexa Fluor® 750, RDye® 750, CF® 750 Dye, DyLight® 750
Excitation/Emission Maximum	758 nm / 775 nm
Shipping Conditions	Ambient temperature
Shipping Instructions	Ambient temperature

ABS/EM SPECTRA

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

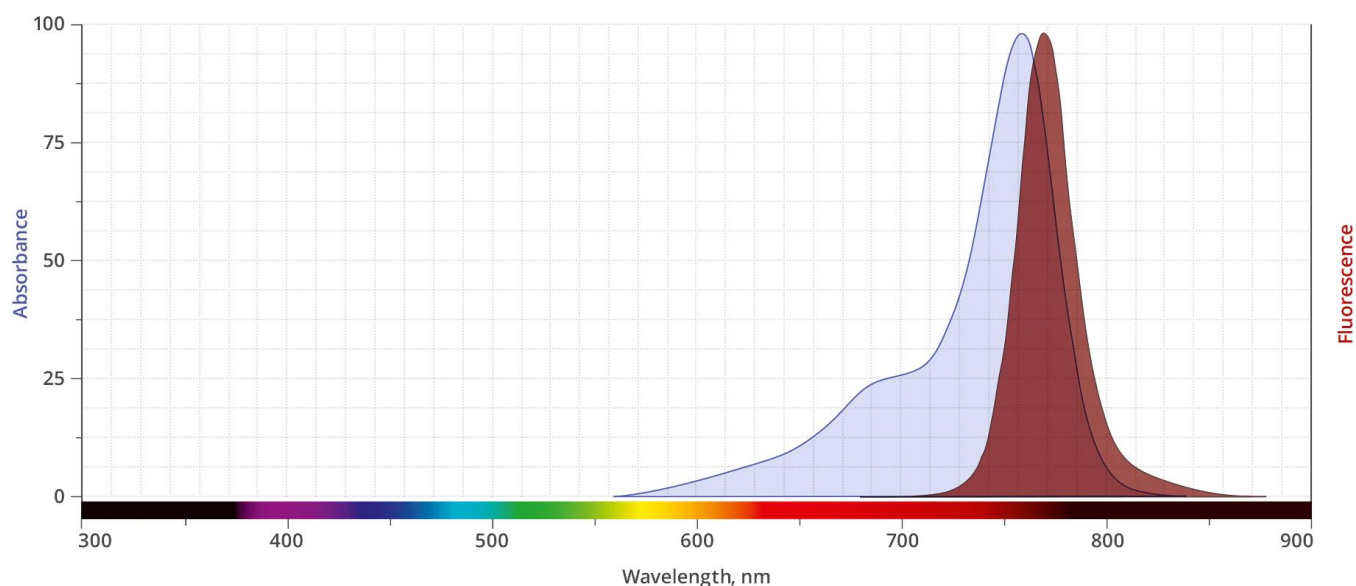

SELECTED REFERENCES

1. Jiang, H., *et al.* (2014). Monitoring Dynamic Glycosylation in Vivo Using Supersensitive Click Chemistry. *Bioconjugate Chem.*, **25**, 698-706. [[PubMed](#)]
2. Uttamapinant, C., *et al.* (2012). Fast, Cell-Compatible Click Chemistry with Copper-Chelating Azides for Biomolecular Labeling. *Angew. Chem. Int. Ed.*, **51**, 5852-56. [[PubMed](#)]
3. Gaebler, A., *et al.* (2016). A highly sensitive protocol for microscopy of alkyne lipids and fluorescently tagged or immunostained proteins. *J. Lipid. Res.*, **57**, 1934-47. [[PubMed](#)]

DOCUMENTS

- [Safety Data Sheet](#)
- [Download CoA](#)
- [Datasheet](#)

GALLERY IMAGES

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

