





ISOTAG BIOTIN AZIDE PACK

SKU: CCT-1501

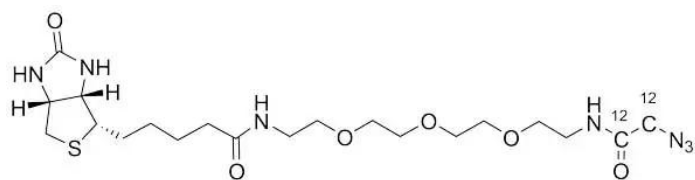

IsoTag L Biotin Azide

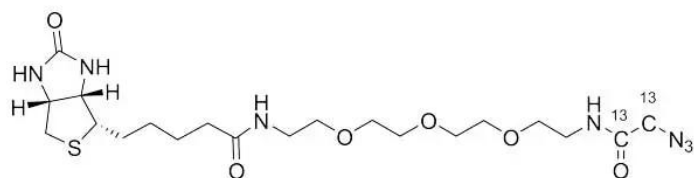

IsoTag H Biotin Azide

DESCRIPTION

IsoTaG Biotin Azide probes take advantage of the mass shift of stable heavy isotope labels (SHLs) to enable mass-independent chemical proteomics platform that helps to address unique challenges of the proteome characterization. In this approach a unique isotopic signature is embedded exclusively into the peptides and it serves as a computationally recognizable full-scan MS reporter. A computational algorithm, termed isotopic signature transfer and mass pattern prediction (IsoStamp), for the detection of recoded species in full-scan mass spectra, was also developed by the Carolyn Bertozzi group. IsoTaG H/L Biotin Azide pack contains 2 mg of light [M] and heavy [M+2] probe encoded by 2 C12 or 2 C13 atoms.

SPECIFICATIONS

Molecular Weight

503.6

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.



Appearance	Oil to grey amorphous solid
Unit Size	1 Pack
Solubility	DMSO, DMF
Storage Instructions	-20C
IsoTag H Biotin Azide	1 vial, 2 mg
IsoTag L Biotin Azide	1 vial, 2 mg
Purchase Note	Add some note
Shelf Life	3 years at -20C
Shipping Instructions	Ambient temperature

SELECTED REFERENCES

1. Woo, C.M., *et al.* (2015). Isotope-targeted glycoproteomics (IsoTaG): a mass-independent platform for intact N- and O-glycopeptide discovery and analysis. *Nat Methods.*, **12**: 561-7. [[PubMed](#)]
2. Woo, C.M., *et al.* (2017). Development of IsoTaG, a Chemical Glycoproteomics Technique for Profiling Intact N- and O-Glycopeptides from Whole Cell Proteomes. *J. Proteome Res.*, **16**: 1706–18. [[PubMed](#)]
3. Goa, G., *et al.* (2017). Small Molecule Interactome Mapping by Photoaffinity Labeling Reveals Binding Site Hotspots for the NSAIDs. *J. Am. Chem. Soc.*, **140**: 4259–68. [[PubMed](#)]
4. Weerapana, E., *et al.* (2010). Quantitative reactivity profiling predicts functional cysteines in proteomes. *Nature* **468**: 790–5. [[PubMed](#)]
5. Woo, C.M., *et al.* (2017). Mapping and Quantification of Over 2000 O-linked Glycopeptides in Activated Human T Cells with Isotope-Targeted Glycoproteomics (Isotag). *Mol. Cell. Proteomics* **17**: 764-75. [[PubMed](#)]

DOCUMENTS

- [Safety Data Sheet](#)
- [Datasheet](#)

GALLERY IMAGES


For research use only. Not intended for therapeutic or diagnostic use in animals or humans.




www.vectorlabs.com

Email: customerservice@vectorlabs.com

Telephone: [\(650\) 697-3600](tel:(650)697-3600)

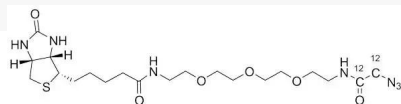

IsoTag L Biotin Azide

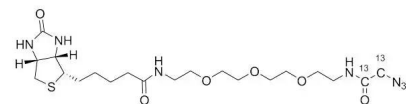

IsoTag H Biotin Azide

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

[IsoTaG Biotin Azide Pack](https://vectorlabs.com/products/isotag-biotin-azide-pack/)

<https://vectorlabs.com/products/isotag-biotin-azide-pack/>