




Email: customerservice@vectorlabs.com

Telephone: [\(650\) 697-3600](tel:6506973600)

SKU: OBD-11643

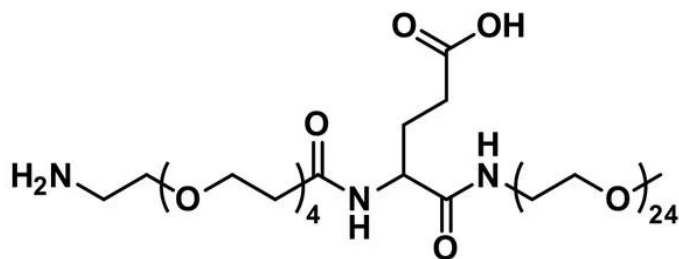

NH2-dPEG®4-Glu(OH)-NH-m-dPEG®24, product number QBD-11643, is one of Vector Laboratories' unique, patented class of modular, designable payload delivery reagents called Sidewinder™. Sidewinder™ products are built on a discrete PEG (dPEG®) backbone for use in antibody-drug conjugates (ADCs) and related constructs. The payload (cytotoxin, dye, small molecule) loads onto the chain's sidearm. The distal end's methoxy-terminated dPEG®24 spacer protects payloads and modifies performance.

The free primary amine of NH₂-dPEG® 4-Glu(OH)-NH-m-dPEG® 24, QBD-11643, reacts with carboxylic acids and their active esters (TFP esters, NHS esters) to form stable amide bonds. It is the attachment point for the Sidewinder™ to the biomolecule (e.g., mAb, Fab, Fab', targeting peptide, and so forth). The backbone contains one carboxylic acid sidearm that reacts with free amines. Following reaction, protection, or functionalization of the attachment point amino group, the carboxylic acid sidearm can be conjugated to an amine-containing payload. The conjugation can be done directly using EDC or another suitable carbodimide. Alternatively, the carboxylic acid can be functionalized with N-hydroxysuccinimide (NHS) or 2,3,5,6-tetrafluorophenol (TFP) to form the active ester, which can then be reacted with a free amine-containing payload.

Sidewinder™ products are designed to facilitate the creation of stable, high-DAR ADCs. Published research has shown that putting a hydrophobic payload close to the antibody surface and protecting it with a SuperHydrophilic™ dPEG® construct is better by many measures of

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

NH₂-dPEG®₄-Glu(OH)-NH-m-dPEG®₂₄

<https://vectorlabs.com/products/nh2-dpeg4-gluoh-nh-m-dpeg24/>



efficacy than putting the payload at the distal end of the linker. This molecule can also modify and optimize BD, cell trafficking and internalization, serum half-life, and immunogenicity.

The dPEG® linkers and spacers in the Sidewinder™ construct are uniform, single molecular-weight PEGs with discrete chain lengths. In contrast, traditional, non-uniform polymer PEG linkers and spacers have a dispersed range of PEG chain lengths, each with a unique molecular weight. Unlike dispersed polymer PEGs, dPEG® products are high-purity compounds with reproducible purity profiles.

Sidewinder™ molecules are also fully designable. NH₂-dPEG®4-Glu(OH)-NH-m-dPEG®24 can be modified to change the spacer lengths, add more sidearm attachment points, add different sidearm attachment points to carry payloads with other reactivities, change the amine attachment group to another reactive group, and many more customizations. Please inquire about your specific needs.

SPECIFICATIONS

CAS Number	N/A
Molecular Weight	1464.72; single compound
Chemical Formula	C ₆₅ H ₁₂₉ N ₃ O ₃₂
Purity	> 95%
Unit Size	50 mg, 250 mg
Solubility	Methylene Chloride, Methanol or Acetonitrile.
Spacers	dPEG® Spacer is 93 atoms and 107.1 Å
Storage Instructions	-20°C; Always let come to room temperature before opening; be careful to limit exposure to moisture and restore under an inert atmosphere; stock solutions can be prepared with dry solvent and kept for several days (freeze when not in use). dPEG® pegylation compounds are generally hygroscopic and should be treated as such. This will be less noticeable with liquids, but the solids will become tacky and difficult to manipulate, if care is not taken to minimize air exposure.
Shipping Instructions	Ambient


For research use only. Not intended for therapeutic or diagnostic use in animals or humans.



DOCUMENTS

- [Safety Data Sheet](#)
- [Datasheet](#)

GALLERY IMAGES

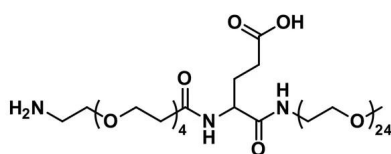

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.