
www.vectorlabs.com
Email: customerservice@vectorlabs.com
Telephone: (650) 697-3600

VECTOR VIP SUBSTRATE KIT, PEROXIDASE (HRP)

SKU: SK-4600

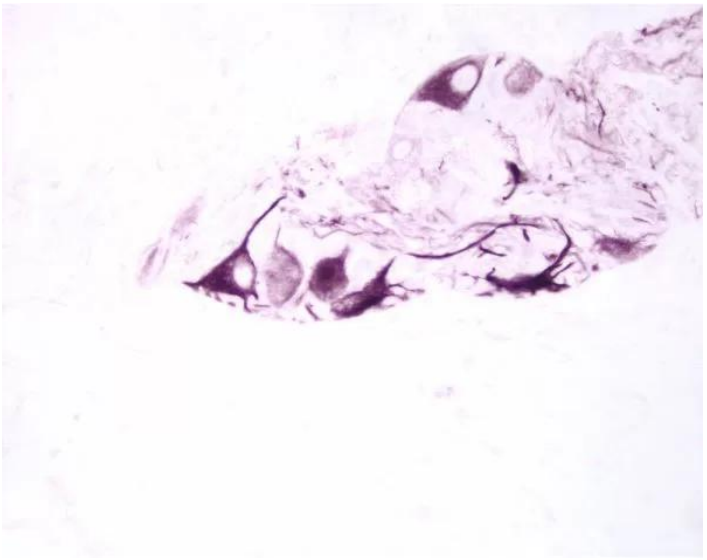

DESCRIPTION

The Vector VIP HRP substrate produces an intense, violet (purple) reaction product, and can be used as an alternative to DAB or as a second color for multiple antigen labeling.

Features:

- Greater sensitivity than conventional substrates
- Consistent and reliable
- Ideal for IHC/ICC/ISH and blots
- Permanent mounting
- Stock solutions supplied in convenient dropper bottles promoting ease of handling
- No wait times for mixing and dissolving powders or tablets
- Sufficient reagents to produce 300 ml of working solution

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

www.vectorlabs.com
Email: customerservice@vectorlabs.com
Telephone: (650) 697-3600

SPECIFICATIONS

Chromogen Color	Purple
Detection Enzyme	Peroxidase
Microscopy	Brightfield, Darkfield, Electron, Spectral Imaging
Mounting Medium Type	Non-Aqueous (Permanent)
Unit Size	1 Kit
Applications	Immunohistochemistry / Immunocytochemistry, In situ hybridization, Blotting Applications, Elispot
Additional Substrate Properties	Contrast with Tissue Pigments, Multiple Labeling
Substrate Reaction Product	Precipitating Reaction Product

KIT COMPONENTS

- 7 ml Vector VIP Reagent 1
- 5 ml Vector VIP Reagent 2
- 5 ml Vector VIP Reagent 3
- 7.2 ml Vector VIP Reagent 4

TECHNICAL INFORMATION


Vector VIP and ImmPACT VIP chromogens also provide excellent color contrast in pigmented tissues such as melanoma or retina (1). Sections stained with either substrate can be viewed by darkfield and electron microscopy. Vector VIP and ImmPACT VIP substrates can be used for both manual and automated staining methods. Stained sections should be dehydrated, cleared, and permanently mounted.

Vector VIP Substrate Kit contains stock solutions in convenient dropper bottles. The sensitivity of Vector VIP substrate is equivalent to Vector DAB Substrate. This product can also be used on blots. The kit provides all of the necessary reagents to prepare about 300 ml of working solution. When stored at 4° C, this kit is stable for one year.

Product review on Biocompare.com.

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

1. Duncan SM, Seigel GM. High-contrast enzymatic immunohistochemistry of pigmented tissues. *J Biol Methods*2016;3(3):e47. doi: 10.14440/jbm.2016.122

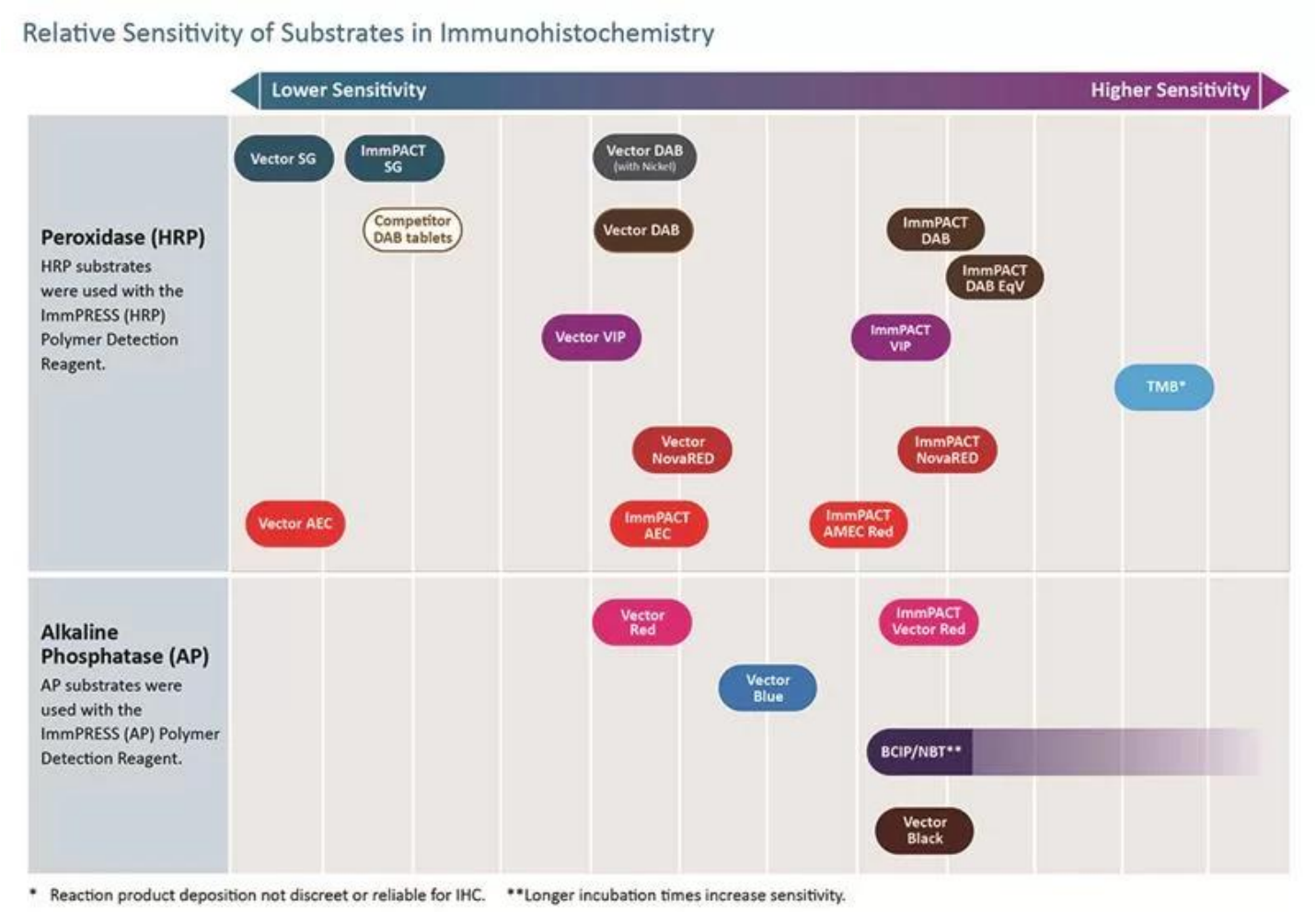

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

www.vectorlabs.com
Email: customerservice@vectorlabs.com
Telephone: (650) 697-3600

Counterstain/Substrate Compatibility Table

This table is designed as a reference to determine the optimal counterstain/substrate combination for your application. Considerations should be given to tissue type, antigen unmasking protocol and other detection parameters to achieve the desired staining intensity.

Substrate	Catalog Number	Vector Hematoxylin and Hematoxylin QS H-3401 and H-3404	Vector Methyl Green H-3402	Vector Nuclear Fast Red H-3403
ImmPACT DAB (brown)	SK-4105	Excellent Contrast	Excellent Contrast	Fair Contrast
ImmPACT DAB EqV	SK-4103	Excellent Contrast	Excellent Contrast	Fair Contrast
DAB (brown)	SK-4100	Excellent Contrast	Excellent Contrast	Fair Contrast
DAB-Ni (gray-black)	SK-4100	Excellent Contrast	Fair Contrast *	Good Contrast
ImmPACT AEC (red)	SK-4205	Excellent Contrast	Counterstain Incompatibility **	Color Incompatibility
ImmPACT AMEC Red (red)	SK-4285	Excellent Contrast	Counterstain Incompatibility **	Color Incompatibility
AEC (red)	SK-4200	Excellent Contrast	Counterstain Incompatibility **	Color Incompatibility
TMB (blue)	SK-4400	Color Incompatibility	Counterstain Incompatibility	Excellent Contrast
ImmPACT VIP (purple)	SK-4605	Fair Contrast	Excellent Contrast	Poor Contrast
Vector VIP (purple)	SK-4600	Fair Contrast	Excellent Contrast	Poor Contrast
ImmPACT SG (blue-gray)	SK-4705	Poor Contrast	Good Contrast	Excellent Contrast
SG (blue-gray)	SK-4700	Poor Contrast	Good Contrast	Excellent Contrast
ImmPACT NovaRED (red)	SK-4805	Excellent Contrast	Excellent Contrast ***	Color Incompatibility
Vector NovaRED (red)	SK-4800	Excellent Contrast	Excellent Contrast ***	Color Incompatibility
ImmPACT Vector Red (magenta)	SK-5105	Excellent Contrast	Excellent Contrast	Color Incompatibility
Vector Red (magenta)	SK-5100	Excellent Contrast	Excellent Contrast	Color Incompatibility
Vector Black (black)	SK-5200	Excellent Contrast	Excellent Contrast *	Excellent Contrast
Vector Blue (blue)	SK-5300	Color Incompatibility	Good Contrast	Excellent Contrast
BCIP/NBT (indigo)	SK-5400	Color Incompatibility	Excellent Contrast *	Excellent Contrast

* This substrate shows a slight decrease in sensitivity following the methyl green protocol. This decrease can be minimized by reducing the heat incubation and acetone rinse times in the methyl green protocol.

** Substrate dissolves in acetone wash.

*** A slight color change in ImmPACT NovaRED™ and Vector NovaRED™ reaction product may be seen using methyl green.

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

www.vectorlabs.com
Email: customerservice@vectorlabs.com
Telephone: (650) 697-3600

Enzyme Substrate Combinations

Recommended combinations of substrates and the recommended order in which they should be used.

First Substrate \ Second Substrate	Alkaline Phosphatase			Peroxidase					
	ImmPACT Vector Red & Vector Red (magenta) SK-5105 & SK-5100	Vector Blue (blue) SK-5300	BCIP/NBT (indigo) SK-5400	ImmPACT VIP & Vector VIP (purple) SK-4605 & SK-4600	ImmPACT DAB, ImmPACT DAB EqV & DAB (brown) SK-4105, SK-4103, SK-4100	DAB-Ni (gray-black) SK-4100	ImmPACT NovaRED & Vector NovaRED (red) SK-4805 & SK-4800	ImmPACT SG & SG (blue-gray) SK-4705 & SK-4700	ImmPACT AEC, ImmPACT AMEC Red & AEC (red) SK-4205, SK-4285, SK-4200
Alkaline Phosphatase									
ImmPACT Vector Red & Vector Red (magenta) SK-5105 & SK-5100		–	–	–	+	+	–	+	–
Vector Blue (blue) SK-5300	+		–	+	+	+	+	+	+
BCIP/NBT (indigo) SK-5400	+	–		+	+	+	+	+	+
Peroxidase									
ImmPACT VIP & Vector VIP (purple) SK-4605, SK-4600	–	+	–		+	+	–	+	–
ImmPACT DAB, ImmPACT DAB EqV & DAB (brown) SK-4105, SK-4103, SK-4100	+	+	+	+		–	–	+	+
DAB-Ni (gray-black) SK-4100	+	–	–	+	+		+	–	–
ImmPACT NovaRED & Vector NovaRED (red) SK-4805, SK-4800	–	+	+	–	+	+		+	–
ImmPACT SG & SG (blue-gray) SK-4705, SK-4700	+	–	–	+	+	–	–		+
ImmPACT AEC, ImmPACT AMEC Red & AEC (red) SK-4205, SK-4285, SK-4200	–	–	–	–	+	–	–	+	

\+ Indicates good contrast – Indicates incompatiblity of substrates for various reasons

CITATIONS

Powered by Bioz © 2023 See more details on Bioz


For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

www.vectorlabs.com
Email: customerservice@vectorlabs.com
Telephone: (650) 697-3600

DOCUMENTS


- Safety Data Sheet
- Enzyme Substrates Brochure
- Enzyme Substrate Combinations Chart
- Counterstain/Substrate Compatibility Chart
- User Guide
- Download CoA
- Datasheet


GALLERY IMAGES

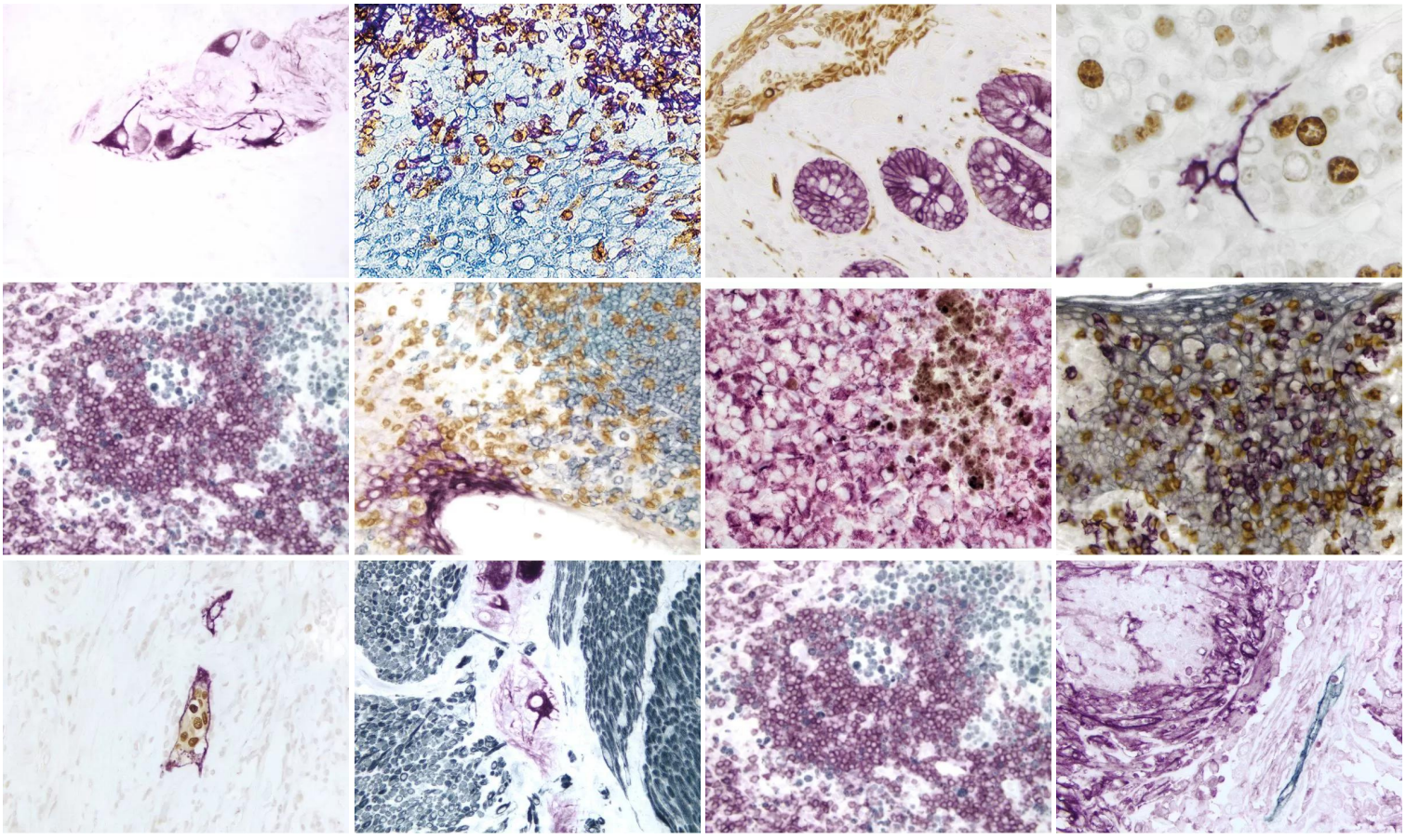



For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

www.vectorlabs.com
Email: customerservice@vectorlabs.com
Telephone: (650) 697-3600

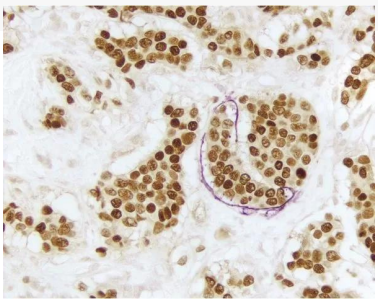

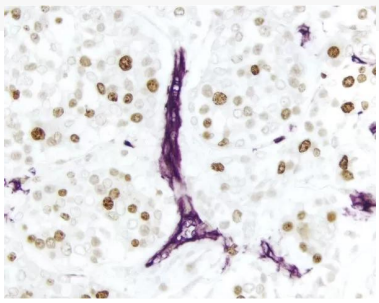

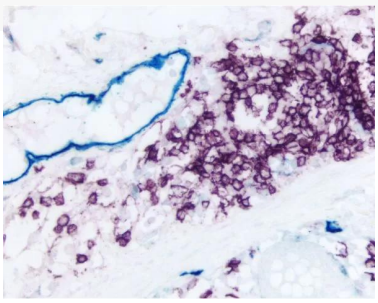

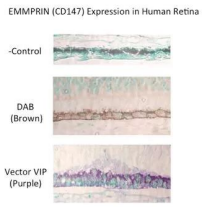



For research use only. Not intended for therapeutic or diagnostic use in animals or humans.